## **Simulation Geometry in GEANT4**

Email: zhaor25@mail2.sysu.edu.cn

School of Physics

- Brief Introduction
- 2 physcis and scoring
- **3** simulation output
- 4 summary

## **Priliminary Geometry Design**

- **1** scintillator cubes  $10 \times 10 \times 10$
- 2 flat film as neutron detector: 4layers
- 3 six light guide arrays
- 4 six PMT [SiPM] arrays





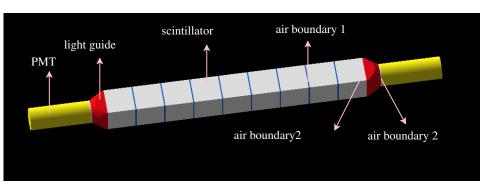


图: scintillator cube

图: scintillator cube+light guide+PMT


# **Details about the Geometry Set-up**

The structure of one dimention detector:



**\bigsilon**: Geometry structure: PMT-lightguide-scintillator-air boundary between [scintillators; scintillator and lightguide; lightguide and PMT cathode]

## **Details about the Geometry Set-up**

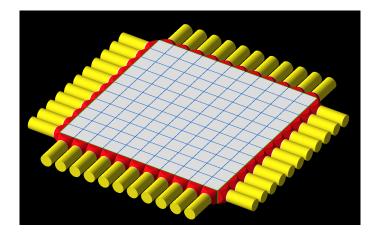

000000000

# **Details about the Geometry Set-up**

Not finished yet.


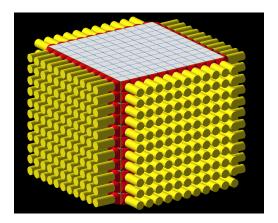


图: 3 dimention detector layout

0000000000


## **Details about the Geometry Set-up**

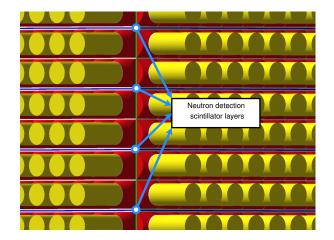



**\bigsilon**: neutron detection scintillator layers in the y direction.


#### flexiable size adjustment

Easy to change the full detector size according to experimental requirements.

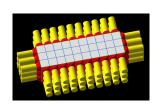

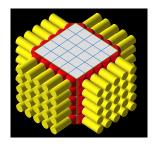

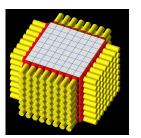

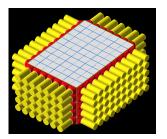


#### THERMAL NEUTRON DETECTOR

- the neutron detector EJ-426.
- flat white thin sheet, 6LiF: (ZnS:Ag)

#### detection princeple:

$$^{6}\text{Li} + ^{1}\text{n} \rightarrow ^{3}\text{H} + ^{4}\text{He} + 4.78\text{MeV}$$
 (1)

The resulting triton and alpha particle are detected by ZnS:Ag phosphor with the broad blue fluorescent spectrum.

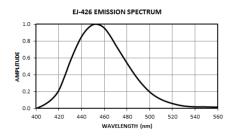

| DETECTION PROPERTIES                             |               |                         |                         |  |
|--|---------------|-------------------------|-------------------------|--|
| Screen Type                                      |               | EJ-426-0                | EJ-426HD2               |  |
| <sup>6</sup> LiF:ZnS Mass Ratio                  |               | 1:3                     | 1:2                     |  |
| <sup>6</sup> Li Density (atoms/cm <sup>3</sup> ) |               | 8.81 × 10 <sup>21</sup> | 1.39 × 10 <sup>22</sup> |  |
| Theoretical<br>N™ Efficiency                     | 0.32 mm thick | 0.23                    | 0.34                    |  |
|  | 0.50 mm thick | 0.34                    | 0.48                    |  |

- choose the formula: EJ-426-0 or EJ-426HD2?
- 2 switch the thickness: 0.32mm or 0.5mm?
- 3 sheet size: 60mm × 60mm?
- 4 do we need backing material?


| BACKING                  |  |        |  |  |
|--------------------------|--|--------|--|--|
| MATERIAL TYPE            | DESCRIPTION                                | SUFFIX |  |  |
| Aluminum Foil            | 50 μm thick foil                           | (none) |  |  |
|                          | 0.25 mm thick sheet                        | -PE    |  |  |
| Clear Polyester Sheet    | Laminated between two 0.25 mm thick sheets | -PE2   |  |  |
| Aluminized Mylar         | 0.12 mm thick sheet                        | -AM    |  |  |
| Pure Aluminum            | 0.5mm thick plate                          | -PA    |  |  |
| High Reflective Aluminum | 0.4mm thick plate                          | -A     |  |  |

#### next to be done

- detector construction.
  - add remain geometry [lightguides and PMTs].
  - attach correct material to each logical volume.
  - other components
- adjustment of physis list
  - about scitillator material and their optical properties
  - optical performance of lightguides
  - response of PMT [SiPM]
  - optical boundaries
- add different primary paticle sources
  - alter the particle type, position, momentum, energy etc.
  - use gps to control theparticle source
- sensitive detector and scoring
- more useractions for output and analyze.

## material of detector components

gamma scintillator:EJ-200Base: Polyvinyl toluene formula: [CH2CH(C6H4CH3)]n

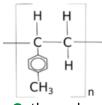

Density: 1.023 g/cm<sup>3</sup> Refraction Index: 1.58

Light Output: No change from -60 $^{\circ}\mathrm{C}$  to 20  $^{\circ}\mathrm{C}$ 

thermal neutron scintillator:EJ-426HD2

<sup>6</sup>LiF: ZnS MassRatio1: 3

 $^6$ Li Density(atoms/cm $^3$ ): 8.81 $\times 10^{21}$ 

<sup>6</sup>Li enriched to minimum of 95 atom percent.

## sensitive detector



- almost finish simple detector geometry.
- other parts of simulation program still in progress.

